

Model 368A

High-Speed, DC Stable Noncontacting Electrostatic Voltmeter



Trek's Model 368A High-Speed Electrostatic Voltmeter provides precision noncontacting measurement of the surface voltages associated with electrostatic imaging processes as employed by the electrophotographic reproduction and laser printing industries. Highly accurate electrostatic voltage measurements are also associated with applications such as photoreceptor research and semiconductor processing. The 368A is used in research and development applications and to measure electrostatic potentials on objects or surfaces such as polymer, film, fabric, and paper.

The miniature sensor probes utilize a patented dual chopper design. This design provides high-speed, DC stable measurements with exceptional low noise and drift performance in wide temperature range conditions, under high humidity, and even in the presence of contaminating particulates. Sensor probe options include standard or high-temperature probes with either an end-viewing or side-viewing orientation.

The voltage of moving or stationary surfaces can be easily measured with high accuracy over a wide range of sensor probe to measured surface spacings. The noncontacting measurement ensures that no charge transfer will take place to modify or destroy the actual data.

A user selectable dual range precision voltage monitor provides a low-voltage replica of the measured electrostatic voltage for external monitoring purposes or for use as a feedback signal in a closed loop control system.

Rack mount options are available to permit up to four independent measurement channels in a standard full rack or for compact bench top operation.

- Noncontacting measurement range
 0 to ±2 kV or peak AC
- Speed of Response less than 200 μs for a 1 kV step
- Accuracy better than 0.1% of full scale
- Easy to read 3½ Digit
 LED display
- Superb noise and drift performance
- Miniature probes
- Precision voltage monitor provides lowvoltage replica of measured voltage
- Up to four independent channels in full rack mountable configuration
- **.** < ∈ compliant



Model 368A Specifications

All specifications are with a probe to surface separation of 2 mm.

Performance

Measurement Range

0 to ±2 kV DC or peak AC.

Measurement Accuracy Voltage Monitor Output

Better than ±0.1% of full scale.

DPM Voltage Display

Better than ±0.15% of full scale, ±1 digit referred to the voltage monitor.

Speed of Response (10% to 90%)

Less than 200 μs for a 1 kV step change.

Stability

Drift with Time

Less than 100 ppm/hour, noncumulative.

Drift with Temperature

Less than 200 ppm/°C.

Probe-to-Surface Separation


2 mm ±1 mm (recommended).

Probe (standard) Aperture Size

3800E-2 end viewing 1.85 mm diameter. 3800S-2 side viewing

2.35 mm diameter. Cable Length

3 m (10 ft), nominal.

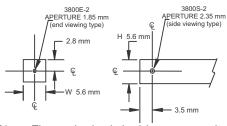

Note: The probe body is driven to match the measured voltage value and must be insulated from ground.

Features

Zero Control

A 10-turn control to null offsets that may be present at the voltage monitor output when the Model 368A is measuring zero volts.

Voltage Display

3 ½ digit LED display.

Range

0 to ±1999 V.

Resolution

1 V.

Zero Offset

Better than or equal to ±1 count, referred to the voltage monitor.

Sampling Rate

2.5 readings per second.

Voltage Monitor Output

A BNC provides a low-voltage replica of the measured voltage.

Scale Factors

1/200th or 1/1000th of the measured voltage (switch selectable). Calibration potentiometers for both scales are accessible from the front panel.

Output Noise

Less than 25 mV rms (measured using the true rms feature of the Hewlett Packard Model 34401A digital multimeter).

Output Impedance

Less than 0.1Ω .

Power ON/OFF

Two position toggle switch.

General

Voltage Monitor Connector

BNC connector.

General (cont.)

Dimensions

Stand-Alone Instrument

133 mm H x 113 mm W x 341 mm D (5.3" H x 4.5" W x 13.5" D).

Half Rack (holds up to two channels) 133 mm H x 271 mm W x 353 mm D (5.3" H x 10.7" W x 13.9" D)

Full Rack (holds up to four channels)
133 mm H x 483 mm W x 353 mm D
(5.3" H x 19" W x 13.9" D).

Weight

Stand-Alone Instrument

2 kg (4.4 lb).

Half Rack

1.5 kg (3.3 lb), add 1.5 kg (3.3 lb) per channel.

Full Rack

2.5 kg (5.5 lb), add 1.5 kg (3.3 lb) per channel.

Ground Receptacle

Banana jack.

AC Line Cord Receptacle

Standard IEC 320 three-prong AC line connector with an integral fuse holder.

Power Requirements Line Voltage

Factory set for one of two ranges: 90 to 127 V AC or 180 to 250 V AC, at 50/60 Hz

(specify when ordering).

Power Consumption

35 VA, maximum.

Operating Conditions Temperature

0°C to 35°C.

Relative Humidity

To 85%, noncondensing. Altitude

To 2000 meters.

Copyright © 2011 TREK, INC 1122/DEC All specifications are subject to change.

Model 368A Ordering Information

IVIOGEI 368A		Probes - Probes	
<u>Item</u>	Part No.*	l <u>ltem</u>	Part No.
Stand-Alone Instrument	368ASU	Standard	
Half-Rack, One Channel		Model 3800E-2 (end viewing)	
Full-Rack, One Channel	368AD 368AT	High-Temperature (up to 60°C) Model 3870ET-2 (end viewing)	
*Add the suffix L for 90 to 127 V AC operation and the suffix H for 180 to 250 V AC operation.	300АQ	Model 3800EC-2 Probe Extension Cable	17385



