

Benchmark

INCU-SHAKER™ Series


The new Incu-Shaker™ series of shaking incubators provide digitally regulated temperature with precisely controlled, circular (orbital) mixing, optimized for the culture of a wide variety of cell types.


The Incu-Shaker Mini is the most compact shaking incubator in its class. Despite the remarkably small footprint, a large 11.5 x 9.5 in. work space accepts several platforms for use with flasks up to 1 liter, tubes and other common vessels. In addition, this model features a ful color, TOUCH screen display for increased control over all parameters

With similar features to the Incu-Shaker Mini, the Incu-Shaker 10L and 10LR models are equipped with a large 18 x 18 inch platform for high capacity applications. (See Chart for a complete listing of the maximum flask capacity.)


The standard platform (included with all three models) features a non-slip, rubber coated surface, ideal for tissue culture flasks, Petri dishes and staining trays. Convenient MAGic Clamp™ platforms, flask clamps and tube racks are also available for use with Erlenmeyer flasks and test tubes. This unique, magnetic attachment method is the easiest way to instantly change between flask clamps of different sizes. Each clamp utilizes a powerful (patent pending) magnetic base to securely attach the clamp to the shaking platform without the use of screws, bolts or tools.

Temperature, shaking speed and run time are displayed on the large LED control panel located on the front of the unit (or Touch screen display on the mini model). The integral microprocessor includes a constant monitoring system that verifies and maintains accuracy through the duration of the program. Sophisticated over-temperature and over-speed ensure long life, safety and sample integrity.

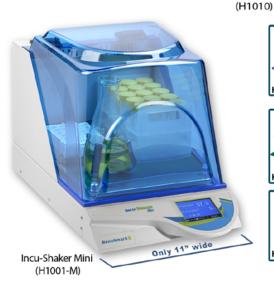

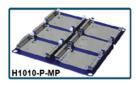



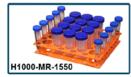

Technical Data:

Temperature Range (Mini): Temperature Range (10L): Temperature Range (10LR): Temperature Uniformity Temperature Accuracy: Shaking Speed Range: Shaking Orbit: Platform Dimensions: Incu-Shaker Mini Incu-Shaker 10L/LR Exterior Dimensions: Incu-Shaker Mini Incu-Shaker 10L Incu-Shaker 10LR

Weight: Incu-Shaker Mini Incu-Shaker 10L Incu-Shaker 10LR Electrical:

Ambient +5 to 70°C (Increments of 0.1°C) Ambient +5 to 60°C (Increments of 0.1°C) Ambient -15 to 60°C (Increments of 0.1°C) +/- 0.25 %

+/-01%

30 to 300rpm (Increments of 1rpm) 19mm (3/4")

9.5 x 11.5 in. / 23.5 x 29.2 cm 17.8 x 17.8 / 45 x 45 cm

11 x 18.8 x 13 in. / 28 x 48 x 33 cm 20 x 26 x 19 in. / 51 x 65 x 47 cm 20 x 31 x 22 in. / 51 X 79 X 56 cm.

*To order an item in 230V, add (-E) to item number.

25 lbs. / 11.5 kg 120 lbs / 55 kg 170lbs / 78 kg 120 or 230V, 50-60Hz



 $n_{\Pi_{\Pi_1}}$



Incu-Shaker 10L

Maximum capacity of Erlenmeyer flasks						
Flask Size:	Incu-Shaker Mini	Incu-Shaker 10L/LR				
50ml	20	49				
125ml	12	36				
250ml	8	25				
500ml	5	16				
1L	2	9				
2L	0	5				

Ordering Information

H1001-M*

Incu-Shaker Mini, includes standard rubber mat platform, 115V Incu-Shaker 10L, includes standard rubber mat platform, 115V Incu-Shaker 10LR, includes standard rubber mat platform, 115V

Optional Platforms and Shelves

	tronno director
H1000-MR	MAGic Clamp universal platform, for Incu-Shaker Mini (9.5x11.5")
H1010-MR	MAGic Clamp universal platform, for Incu-Shaker 10L/LR (18x18")
H1000-P-MP	Microplate platform, 4 position for Incu-Shaker Mini
H1010-P-MP	Microplate platform, 6 position for Incu-Shaker 10L/LR
H1000-P-SP	Spring rack platform, holds uncommon vessels for Incu-Shaker Mini
H1010-P-SP	Spring rack platform, holds uncommon vessels for Incu-Shaker 10L/LR
H1010-SH	Optional Shelf, for Petri Dishes, microplates, trays, etc. (14 x 5.5")
	for Incu-Shaker 10L/LR (max. 3)

Flask & Tube Rack Accessories for MAGic Clamp Platforms Flask Clamp, 50ml (H1000-MR: max. 20 / H1010-MR: max. 49)

		H1000-MR-125	Flask Clamp, 125ml ((H100-MR: max. 12 / H10	10-MR: max. 36)		
		H1000-MR-250	Flask Clamp, 250ml	H1000-MR: max. 8 / H10	10-MR: max. 25)		
		H1000-MR-500	Flask Clamp, 500ml	H1000-MR: max. 5 / H10	10-MR: max. 16)		
		H1000-MR-1000	Flask Clamp, 1000ml	(H1010-MR: max. 9)			
		H1000-MR-2000	Flask Clamp, 2000ml	(H1010-MR: max. 5)			
		H1000-MR-MP	Microplate Clamp, (H1	1000-MR: max. 4 / H1010	-MR: max. 6)		
		H1000-MR-1550	Tube Rack, 20x50 & 30	x15 ml tubes (H1000-MR: m	nax.1 / H1010-MR: max. 4)		
		H1000-MR-T15	Tube Rack, 32x15 tube	es (H1000-MR: max. 2 / H1	010-MR: max. 8)		
		H1000-MR-T50 Tube Rack, 15x50 tubes (H1000-MR: max. 2 / H1010-MR: max. 8) H1000-MR-T600 Bottle Rack, 4x500/600ml conical bottles (H1000-MR: max. 1 / H1010-MR: max. 2)					
	H1000-MR-TSR MAGic Clamp™ tilted holder for 15mL and 50mL tube racks						
	lb.	H1000-MR-CMB	Adjustable flask clamp	for MAGic clamp platform,	accepts 125ml (max. 36),		
	IIIIII II ii .		250ml (max. 25) or 500	0ml (max. 16)			
	241111			Additional accessories available. Pleas	e visit www.BenchmarkScientific.com for more details		
		HIII.		Took Familians and	OO Weekington Charact	Ī	
		SHIIIII III III III III III III III III		Test Equipment	99 Washington Street		
		- 1111111111111111111111111111111111111	r.	Depot	Melrose, MA 02176		
		.,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,	H1b.	- VIV	Phone 781-665-1400		
			HIIIII Do.	1-800-517-8431	Toll Free 1-800-517-8431		
ice	€.		944444444	1-000-011-0401	10111100 1 000 317 0 131		
			.,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,,	A vertice at the second	Took Facilians and Daniel and		
					.TestEquipmentDepot.com		
					Little .		

